

MA 1

Super-Lightweight Hearing Screener

The MA1 is MAICO's smallest solution for fast hearing screening. Its sleek design and intuitive handling make testing quick and comfortable. The MA1 is always within reach – anytime, anywhere.

Features at a Glance

- Handheld air-conduction screener
- Easy one hand operation
- Exceptionally lightweight
- Test frequencies: 500, 1000, 2000 and 4000 Hz
- Pure tone signal from 15 dB HL to 50 dB HL
- Battery operated (2 AA batteries)
- Automatic power off
- AC headphones DD45
- Optional AC noise-enclosure headphones DD65 v2



Applications

Designed for pediatricians, childhood screening, and hearing aid dispensers, the MA1 is a reliable handheld screener that enables quick and accurate hearing assessments. Its lightweight, portable design simplifies daily screenings in any setting, supporting early detection and effective hearing care.

Small and Handy

With its super-lightweight design the MA1 is small enough to fit in the palm of your hand. MA1 is our lightest and most handy air conduction hearing screener. Frequencies and intensities can be presented by one hand operation.

Ideal for Fast Screening Purposes

The device is capable delivering a pure tone signal to one ear at a time. The stimulus range is 15 to 50 dB HL in level steps of 5 dB. Test Frequencies range from 500 to 4000 Hz. The device is equipped with the AC headphones DD45.

If the MA1 is used in noisy environment, we recommend to choose the optional AC headphones DD65 v2 with high attenuation.

Extended Battery Life

The MA 1 is powered by 2 AA alkaline batteries. The device turns off automatically when inactive and turns on again when pressing any of the 7 buttons.


MA 1 – Components & Accessories



MA 1 handheld device



AC headphones DD45

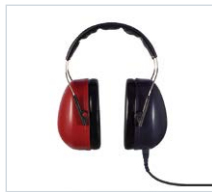

AC headphones DD65 v2




Carrying case



Audiogram pad

Components and accessories vary by region and configuration.

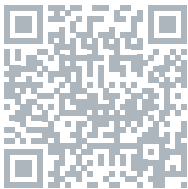

Watch our MA 1 training video

Find out more about other MAICO solutions

easyTone



Fast and Easy
Tablet Audiometry

ERO•SCAN



Easy-To-Use
OAE Handheld Device



Global Contact:

MAICO Diagnostics GmbH

Sickingenstr. 70 - 71 · 10553 Berlin · Germany

Tel.: +49 (30) 70 71 46-50

sales@maico.biz · www.maico-diagnostics.com

US Contact:

MAICO Diagnostics US

10393 West 70th Street · Eden Prairie · MN 55344 · USA

Tel.: +1 (888) 941-4201 · Fax: +1 (952) 903-4100

info@maico-diagnostics.com · www.maico-diagnostics.com



09/2025