

Hearing Screening for All Ages



ARRIVING
SPRING
2027

easyTone **Connect**

Connected hearing screening

easyTone Connect

PC-based hearing screening that brings testing, records and reporting together in one workflow

easyTone Connect is a PC-based solution that combines screening, data management and reporting in one streamlined workstation.

Operated through MAICO Sessions software, it combines full-function pure tone screening with RASA™ to deliver accurate, efficient testing for a connected workflow.

Ideal for school screening programs, occupational health, community and primary care settings.

Powered by the easyTone IEC/ANSI Type 4 audiometer



Founded to Make Hearing Screening Easier

MAICO Diagnostics, founded in 1937 by Leland A. Watson, was established to meet a growing need among doctors for a reliable, scientific tool to assess hearing loss.

From the beginning, our mission has been to create accurate, easy-to-use hearing test equipment that supports early detection of hearing loss and improves quality of life around the world.



Applications

Adult Screening

Designed for occupational health providers, clinical environments and hearing care professionals, easyTone Connect streamlines screening with standardized workflows, digital records and efficient reporting.

Access to audiograms and patient records supports informed care decisions, appropriate referrals and ongoing patient management.

The quick and dependable screening process helps identify potential hearing loss and guides individuals toward the next step in their hearing care journey.

Childhood Screening

Built for large-scale childhood screening, easyTone Connect helps users efficiently assess large groups of children with fast, accurate, and repeatable results.

Built-in noise monitoring helps maintain reliable testing in challenging environments, while standardized workflows, digital records, and streamlined reporting make it easier to manage screening across schools, public health programs, and community initiatives.

Its PC-based operation also provides simple integration into existing IT environments.

Key Features

- PC-based operation for centralized hearing screening
- Rapid Assisted Screening Audiometry (RASA™)
- Pure tone audiometry (manual and automatic)
- Data management and reporting with MAICO Sessions software
- Seamless integration to third-party EMR systems
- Classification of hearing loss
- Room noise monitoring
- Audiogram comparison
- Lightweight and portable



RASA™: Fast, Consistent Hearing Screening

Supported by RASA™ (Rapid Assisted Screening Audiometry), easyTone Connect streamlines mass hearing screening with a standardized, semi-automatic workflow.

Select your Workflow

RASA™ makes it simple to create a customized workflow around your individual screening needs. The built-in protocol wizard guides users step by step through selecting frequencies, Pass/Refer criteria and test sequences.

With only the essential controls and information displayed, users can follow a clear, consistent workflow without unnecessary complexity.

Easy Interpretation

Clear on-screen Pass/Refer results simplify interpretation, helping to reduce errors and support consistent screening.

With a typical screening completed in under two minutes, RASA™ helps users screen more patients quickly and efficiently.

Guided Screening

Designed to save time and simplify hearing screening, RASA™ guides users through a streamlined testing workflow.

Users simply present tones and record results, while the device automatically manages frequency progression and ear switching.


This reduces complexity, minimizes operator input, and helps deliver fast and consistent results.

Confident Results

RASA™ helps users maintain confidence in their test results by monitoring ambient noise during screening.

The headphone's integrated noise monitor identifies background noise that could affect test accuracy, with noise levels displayed via an on-screen indicator.

Screeners can instantly verify testing conditions, supporting accurate results across different environments.

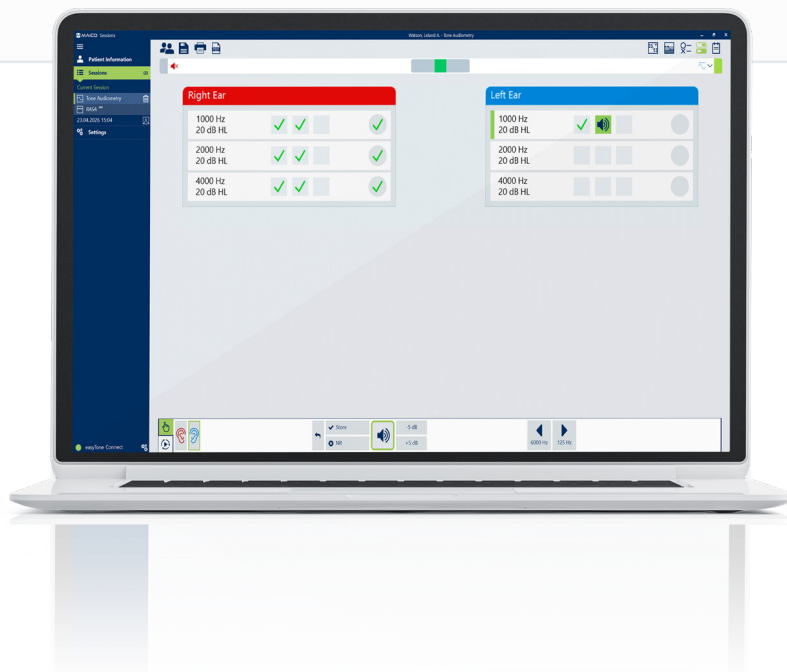

Features & Benefits

Connected Data Management

One workstation - everything connected. easyTone Connect brings screening, patient records and reporting systems together in a single interface.

Results are automatically transferred to integrated databases and EMR systems. These include MAICO Sessions Database, Noah 4, OtoAccess® and XML/GDT-compatible systems, eliminating manual data transfer and duplicate entry. From testing to documentation to referral, everything stays connected.

PC-based

easyTone Connect is operated through MAICO Sessions on a Windows PC.* Simply connect the calibrated headphone to your personal device via USB and start screening.

Patient records, audiograms and reporting tools are all accessible from a single interface, eliminating manual data entry and the need to switch between applications.

Whether used at a dedicated screening station or across multiple locations, easyTone Connect adapts seamlessly to your workflow.

***PC must meet specified system requirements.**

Compact by Design

Unlike traditional audiometry setups, easyTone Connect does not require additional hardware, using only the easyTone headphone audiometer and your existing PC or laptop.

For clinics with limited space or individual workstations, the device fits into existing setups without clutter or complexity.


This means less equipment to store, maintain or transport between rooms.



Trusted easyTone Technology

easyTone Connect uses the same high-attenuation circumaural headphones as other easyTone solutions. Its unique design houses both an audiometer and its calibration, building on proven performance to deliver the next generation of hearing screening.

Components, Accessories & Software

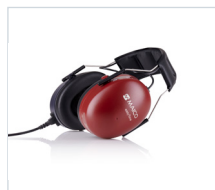

easyTone
audiometer
headphones



Patient response
switch



Carrying case



MAICO Sessions

Components and accessories vary by region and configuration.

Find out more about other MAICO solutions

ERO•SCAN

Easy-to-Use OAE
Handheld Device



touchTymp

Intuitive Touch Screen Tympanometry
with Advanced Testing



Global Contact:

MAICO Diagnostics GmbH

Sickingenstr. 70 -71 · 10553 Berlin · Germany

Tel.: +49 (30) 70 71 46-50

sales@maico.biz · www.maico-diagnostics.com

US Contact:

MAICO Diagnostics US

10393 West 70th Street · Eden Prairie · MN 55344 · USA

Tel.: +1 (888) 941-4201 · Fax: +1 (952) 903-4100

info@maico-diagnostics.com · www.maico-diagnostics.com



09/2026